



OBio...

Supporting Care to Cure

A Division of Orchard Medical



OBio ESR 10-Pro

Green Solution, Reduced Waste, Lower Expenses

**NOT REQUIRED
CITRATE TUBE**




- ESR test performed on same tubes in which blood samples taken (EDTA/CBC)
- Test can be done on small volume of whole blood (1.5mL)

ESR test results on a **compact 7" touch screen**,
designed for seamless interaction and clear visualization



- User friendly & Environment friendly
- Keep on loading of samples & result in 20 minutes
- No waste material, bringing down the cost of liquid waste disposal

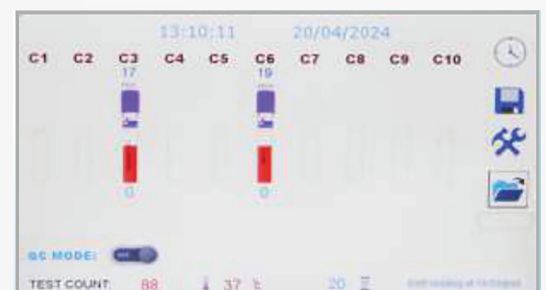

**RESULT CALCULATION
BY USING MODIFIED WESTERGREN METHOD**

OBio ESR 10-Pro

Technical Specifications

Power supply	110-230V AC 50Hz/60Hz
Dimension	189 mm(Length) x 282 mm(Width) x 231 mm(Height)
Weight	4 Kg
Working Temperature	15°C-30°C
Relative Humidity	20 to 80% without condensation
Storage conditions	10°C-40°C
LIS interface	Via RS-232, HL7 protocol
USB Port	2 USB port for printer & barcode reader

Product code

Instrument	OBIO-ESR-I-10	OBIO-ESR 10-Pro
------------	---------------	-----------------

Consumables	OBIO-ESR-R-500	RECHARGE TEST 500 Nos
-------------	----------------	-----------------------

Accessories	OBIO-ESR-BR	OBIO-ESR-TP	OBIO-ESR-TP-P
	BARCODE READER	THERMAL PRINTER	THERMAL PAPER



Orchard Medical
154, Teghoria Main Road Kolkata -700157
Mobile no: - +91 84200 69980
Website- www.orchardmedical.co.in
Email: sales.orchardmedical@gmail.com



OBio...
Supporting Care to Cure
A Division of Orchard Medical